

From a laboratory report to a clearer conversation.

A page-by-page guide to asking a question, reviewing extracted values and exploring the explanation.


This guide captures the real redesigned interface. Mocked stages are marked on the screenshots and in their captions. It is not a record of live provider or clinical validation.



Step 01

Enter the workspace

Open ResultScope. Choose Open workspace or scroll to reach the conversation on the same page. Try an example loads synthetic text without submitting it. The desktop transition follows your normal scroll; reduced-motion preferences and mobile layouts use a static transition.

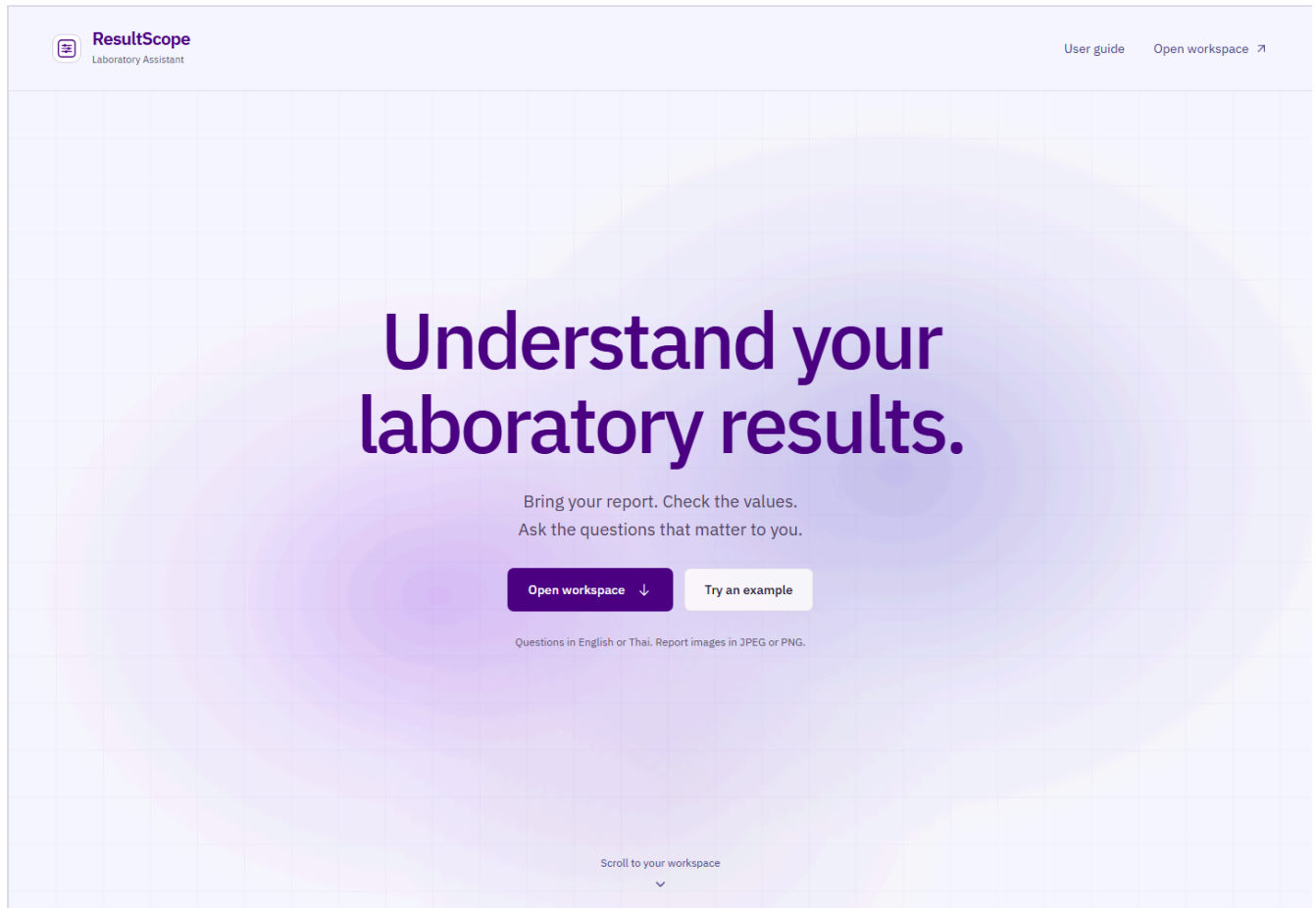

Actual landing page; no provider call.

Step 02

Start in the workspace

In the conversation, select Your result or question. You do not need an administrator account to ask a laboratory question. Type in English or Thai. Include the value, unit and the reference range printed on your report.

The screenshot shows a web interface titled "Laboratory conversation" in the top left corner. In the top right corner, there is a "Report" button with a document icon. Below the title bar, a yellow banner contains the text: "Synthetic business data for demonstration only; no real service is provided." The main area has a light purple background with a grid pattern. In the center, there is a circular icon with a document and a list. Below this icon, the heading "What would you like to understand?" is displayed, followed by the instruction "Ask about a result, or attach your laboratory report." There are two input fields: "Understanding haemoglobin" and "Creatinine and eGFR", each with a right-pointing arrow icon. Below these fields is a larger white box titled "Your result or question" with the placeholder text "Type a question, or paste values with their units and reference ranges...". At the bottom left of this box is an "Attach report" button with a paperclip icon. At the bottom right is a "Send" button with a right-pointing arrow, and a character count "0 / 12,000" is positioned just above it. At the bottom of the interface, there are three lines of small text: "Use the reference range printed on your report." with a "Load example" link; "Examples use synthetic values. No diagnosis or treatment advice."; and "Educational information. No diagnosis, prescriptions or treatment changes." with a "View calculation rules" link.

Real application page; captured with provider network disabled.

Step 03

Try a synthetic example

Select Load example to fill the input. Review the text before choosing Send. Loading an example does not call a provider; submitting may do so only when live access has been explicitly enabled by an operator.

Laboratory conversation

Report

Synthetic business data for demonstration only; no real service is provided.

What would you like to understand?

Ask about a result, or attach your laboratory report.

Understanding haemoglobin

Creatinine and eGFR

Your result or question

CBC: Hb 10.8 g/dL (12-16), MCV 72 fL (80-100), RDW 17.2% (11.5-14.5), Ferritin 7 ng/mL (15-150). Please summarize key findings, relationships, and missing context.

Attach report

163 / 12,000

Send

Use the reference range printed on your report.

Load example

Examples use synthetic values. No diagnosis or treatment advice.

Educational information. No diagnosis, prescriptions or treatment changes.

View calculation rules

Real example control. Example values and ranges are synthetic.

Step 04

Upload and review a report

Choose Attach report or Report to open the document drawer. Select a de-identified JPEG or PNG up to 3 MB. The drawer shows the selected filename and an image preview; reading an image requires an authorized OCR provider. Check each extracted value, unit and range against the image. Edit incorrect fields. If extraction is unavailable, type the values instead. PDF upload is not supported.

Documentation capture · Mock OCR response · Synthetic report · No provider call

Laboratory conversation

Your report

Synthetic business data for demonstration only; no real service is provided.

What would you like to understand?

Ask about a result, or attach your laboratory report.

Understanding haemoglobin

Creatinine and eGFR

Your result or question

Type a question, or paste values with their units and reference ranges...

Attach report

0 / 12,000

Send →

Use the reference range printed on your report.

Load example

Examples use synthetic values. No diagnosis or treatment advice.

Educational information. No diagnosis, prescriptions or treatment changes.

View calculation rules

Documentation capture · Mock OCR response · Synthetic report · No provider call

Your report

Check the image and extracted values.

Upload a laboratory report

JPEG or PNG, up to 3 MB. Remove personal details first. PDF is not supported.

Choose File

No file chosen

synthetic-report.png

34.9 KB · PNG

Review and edit the values

SYNTHETIC LABORATORY REPORT

Documentation fixture only. No patient information.

Test	Value	Unit	Reference range
Haemoglobin	10.8	g/dL	12 - 16
MCV	72	fL	80 - 100
Ferritin	7	ng/mL	15 - 150

Ranges are invented for demonstration, not clinical guidance.

Review values extracted from the image

Text from the image is unverified data. Edit unclear values before confirming.

Hb

Value

Unit

10.8

g/dL

Reference range

12-16

MCV

Value

Unit

72

fL

Reference range

80-100

Real review UI with a MOCK OCR response. No OCR API call was made.

Confirm the fields

Choose Confirm values only after checking every field. The review panel changes to Confirmed extracted values. The confirmed fields become read-only and the question box is filled with a request to explain them. On desktop choose Return to conversation; on mobile the panel closes after confirmation. Choose Send when ready. Before sending, Discard image clears the pending attachment. After sending, start a new analysis to replace a report.

Documentation capture · Mock OCR response · Synthetic report · No provider call

Laboratory conversation

Your report

Synthetic business data for demonstration only; no real service is provided.

What would you like to understand?

Ask about a result, or attach your laboratory report.

Understanding haemoglobin

Creatinine and eGFR

Your result or question

Please explain the values extracted from this laboratory report.

Attach report 64 / 12,000 Send

Use the reference range printed on your report.

Load example

Examples use synthetic values. No diagnosis or treatment advice.

Educational information. No diagnosis, prescriptions or treatment changes.

View calculation rules

Documentation capture · Mock OCR response · Synthetic report · No provider call

Your report

Check the image and extracted values.

Test	Value	Unit	Reference range
Haemoglobin	10.8	g/dL	12 - 16
MCV	72	fL	80 - 100
Ferritin	7	ng/mL	15 - 150

Ranges are invented for demonstration, not clinical guidance.

Confirmed extracted values

These values are confirmed for this analysis. Start a new analysis to use a different report.

Hb

Value

Unit

10.8

g/dL

Reference range

12-16

MCV

Value

Unit

72

fL

Reference range

80-100

Ferritin

Value

Unit

7

ng/mL

Reference range

15-150

Discard image

Return to conversation

Real confirmation controls with a MOCK confirmation endpoint. This capture does not validate extraction accuracy.

Read the result summary

Read the extracted values before the explanation. Select a value to inspect its supplied reference range. Below range, Above range, Within range and Range unknown are comparison states, not diagnoses. Read the explanation together with missing context.

Documentation capture - Fixed mock answer - Synthetic values - No live LLM validation

Laboratory conversation

Your report

Synthetic business data for demonstration only; no real service is provided.

Your conversation

Start a new analysis

Result summary

1 Read 2 Check 3 Explain

Message sent to the system

Show your message

3 numeric values found

Rules 2026.08

Hb
10.8 g/dL
Below range

MCV
72 fL
Below range

Ferritin
7 ng/mL
Below range

Hb
10.8 g/dL

Below range

12 g/dL 16 g/dL

Calculated against the reference range supplied in this message.

Report reference range 12-16 g/dL

Your results in context

Explanation complete

Reading the supplied values

The synthetic example contains haemoglobin **10.8 g/dL** against a supplied range of **12-16 g/dL**. This value is below that supplied range.

The MCV and ferritin values are also below the example ranges. These comparisons do not establish a diagnosis.

What to check next

Confirm the values, units and ranges against the original report. Your laboratory or clinician can explain how your circumstances affect interpretation.

Documentation demonstration only: this explanation is a fixed mock response, not a live AI answer.

Sources used

View source details

Calculation details

Inspect the supplied ranges and calculation rules

Ask a follow-up

Ask about a value, the source, or missing context

☐ View report

Send

For a different report, start a new analysis.

Educational information. No diagnosis, prescriptions or treatment changes.

View calculation rules

Inspect the source

Expand Sources used when it is available. Check the document title, organization, version and page before opening its link. Calculation details explains the deterministic rules. A displayed source is not proof of clinical applicability to your case.

Documentation capture · Fixed mock answer · Synthetic values · No live LLM validation

Laboratory conversation

Your report

Synthetic business data for demonstration only; no real service is provided.

Your conversation

Start a new analysis

Result summary

1 Read 2 Check 3 Explain

Message sent to the systemShow your message

3 numeric values foundRules 2026.08

Hb

10.8 g/dL

Below range

MCV

72 fL

Below range

Ferritin

7 ng/mL

Below range

Hb

10.8 g/dL

12 g/dL16 g/dL

Below range

Calculated against the reference range supplied in this message.

Report reference range12–16 g/dL

Your results in contextExplanation complete

Reading the supplied values

The synthetic example contains haemoglobin **10.8 g/dL** against a supplied range of **12–16 g/dL**. This value is below that supplied range.

The MCV and ferritin values are also below the example ranges. These comparisons do not establish a diagnosis.

What to check next

Confirm the values, units and ranges against the original report. Your laboratory or clinician can explain how your circumstances affect interpretation.

Documentation demonstration only: this explanation is a fixed mock response, not a live AI answer.

Sources usedView source details

Synthetic report – documentation fixture

Demonstration data only

Page 1 - synthetic

License: Project-authored synthetic example; not a clinical source

Open source

Calculation details

Inspect the supplied ranges and calculation rules

Ask a follow-up

Ask about a value, the source, or missing context

View reportSend

For a different report, start a new analysis.

Educational information. No diagnosis, prescriptions or treatment changes.

View calculation rules

Step 08

Ask a follow-up or start again

Use Ask a follow-up to ask about the same results. Choose Start a new analysis before using a different report. Wait for the reset to succeed. Do not mix reports from different people in one conversation.

Documentation capture · First mock answer · Synthetic values · No live LLM validation

Laboratory conversationYour report

Synthetic business data for demonstration only; no real service is provided.

Your conversationStart a new analysis

Result summary1 Read 2 Check 3 Explain

Message sent to the systemShow your message

3 numeric values foundRisks: 2024.08

Hb10.8 g/dLBelow range

MCV72 fLBelow range

Ferritin7 ng/mLBelow range

Hb10.8 g/dL

12 g/dL54 g/dL

Calculated against the reference range supplied in this message.

Report reference range12-16 g/dL

Your results in contextExplanation complete

Reading the supplied values

The synthetic example contains haemoglobin 10.8 g/dL against a supplied range of 12-16 g/dL. This value is below that supplied range.

The MCV and ferritin values are also below the example ranges. These comparisons do not establish a diagnosis.

What to check next

Confirm the values, units and ranges against the original report. Your laboratory or clinician can explain how your circumstances affect interpretation.

Documentation demonstration only: this explanation is a fixed mock response, not a live AI answer.

Sources usedView source details

Synthetic report – documentation fixture

Demonstration data only

Page 1: synthetic

License: Project-authored synthetic example; not a clinical source

Open source

Calculation details

Inspect the supplied ranges and calculation rules

Follow-up in the same context1 Read 2 Check 3 Explain

Message sent to the systemShow your message

3 numeric values foundRisks: 2024.08

Hb10.8 g/dLBelow range

MCV72 fLBelow range

Ferritin7 ng/mLBelow range

Hb10.8 g/dL

12 g/dL54 g/dL

Calculated against the reference range supplied in this message.

Report reference range12-16 g/dL

Your results in contextExplanation complete

Keep the same context

This follow-up uses the synthetic report above. The supplied reference range is the basis of the comparison. Ask the laboratory which range applies to your report.

Documentation demonstration only.

Sources usedView source details

Calculation details

Inspect the supplied ranges and calculation rules

Ask a follow-up

Ask about a value, the source, or missing context

View reportSend

For a different report, start a new analysis.

Education information. No diagnosis, prescriptions or treatment changes.View calculation rules

Real follow-up controls with MOCK SSE payloads. Session behavior is additionally covered by existing offline tests.

Recover when a provider is unavailable

Read the error and keep the values you entered. Provider access may be disabled, unconfigured or out of budget. Ask the operator to resolve the specific condition. Do not repeatedly retry an exhausted budget. Never paste API keys into the question box.

Laboratory conversation

Report

Synthetic business data for demonstration only; no real service is provided.

Your conversation

Start a new analysis

Result summary

1 Read 2 Check 3 Explain

Message sent to the system

Show your message

4 numeric values found

Rules 2026.08

Hb

10.8 g/dL

Below range

MCV

72 fL

Below range

RDW

17.2 %

Above range

Ferritin

7 ng/mL

Below range

Hb

10.8 g/dL

Below range

12 g/dL 16 g/dL

Calculated against the reference range supplied in this message.

Report reference range 12-16 g/dL

The explanation could not be completed

Explanation unavailable

Your entered values remain visible. Read the message below for the next step.

Calculation details

Inspect the supplied ranges and calculation rules

The explanation service has not been configured. Ask the administrator to complete setup.

Try again

Ask a follow-up

Ask about a value, the source, or missing context

View report

Send

For a different report, start a new analysis.

Educational information. No diagnosis, prescriptions or treatment changes.

View calculation rules

Actual offline/missing-configuration route; no browser response mock.

Recognize the scope boundary

ResultScope handles laboratory information. An unrelated question is refused locally. Ask a laboratory question instead. The same boundary excludes diagnosis, prescribing and treatment changes.

Laboratory conversation

Report

Synthetic business data for demonstration only; no real service is provided.

Your conversation

Start a new analysis

Outside laboratory scope

Message sent to the system

Show your message

Rule-based response

Explanation complete

ResultScope is focused on laboratory results.

Share a result with its unit and reference range, such as Hb 10.8 g/dL (12–16), or ask what a result such as HbA1c 6.1% may mean.

Explain this CBC panel

How should I read Creatinine and eGFR together?

What should I read alongside high AST/ALT?

What is HbA1c?

Ask a follow-up

Ask about a value, the source, or missing context

View report

Send →

For a different report, start a new analysis.

Educational information. No diagnosis, prescriptions or treatment changes.


View calculation rules

Actual deterministic scope refusal; no LLM call.

Step 11

Administrator sign in

Operators can open Administrator from the footer. The login is available only in explicit local-demo mode. The demonstration credentials admin / 1234 must never be used online. Ordinary users do not need this account.

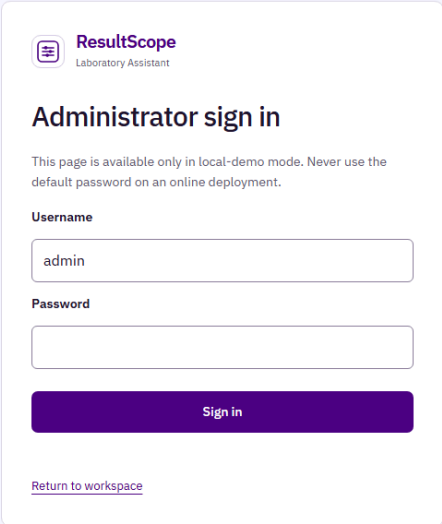

The screenshot shows a web interface for 'ResultScope Laboratory Assistant'. The title is 'Administrator sign in'. Below the title, a warning message states: 'This page is available only in local-demo mode. Never use the default password on an online deployment.' There are two input fields: 'Username' with the value 'admin' and 'Password' which is empty. A blue 'Sign in' button is below the password field. At the bottom, there is a link 'Return to workspace'.

ResultScope
Laboratory Assistant

Administrator sign in

This page is available only in local-demo mode. Never use the default password on an online deployment.

Username

Password

Sign in

[Return to workspace](#)

Actual local-demo login page. The captured environment contains no private credentials.

Step 12

Configure providers privately

Choose the provider and model, then enter a replacement API key only in the password field. Leave it blank to retain an existing key. Save settings stores the configuration without a provider call. Test with mock does not verify live connectivity. SystemOne remains shadow-only.

 **ResultScope**
Laboratory Assistant

Log out

Provider settings

Settings are stored only for local demo use. Existing API keys are never displayed.

OpenTyphoon LLM

OpenAI-compatible chat completions. Endpoint: <https://api.opentyphoon.ai/v1>. Key: not configured.

Provider

Model

Timeout (seconds)

OpenTyphoon LLM (typhoon_llm)

typhoon-v2.5-30b-a3b-instruct

60

Replace key

Leave blank to keep the current key

☒ enabled

Test with mock

Delete key

Typhoon OCR

Separate OCR document contract using multimodal chat payloads. Endpoint: <https://api.opentyphoon.ai/v1>. Key: not configured.

Provider

Model

Timeout (seconds)

Typhoon OCR (typhoon_ocr)

typhoon-ocr

60

Replace key

Leave blank to keep the current key

☒ enabled

Test with mock

Delete key

OpenThai-SystemOne

Typed decision API; no generated text; shadow-only by default. Endpoint: <https://api.iapp.co.th/v3/store/openthai/systemone>. Key: not configured. shadow-only

Provider

Model

Timeout (seconds)

OpenThai-SystemOne (openthai_syste

openthai-systemone

60

Replace key

Leave blank to keep the current key

☐ enabled

Test with mock

Delete key

Save settings

 Save does not call a provider. Tests are separate and use a mock by default.

[Return to workspace](#)

Actual settings page with unconfigured slots. No Save, provider test or live request was performed during capture.

Step 13

Use the mobile workspace

On a phone, choose Open workspace or scroll down from the landing page. Type in the question box and select Send. Report opens a full-screen document panel. Use Close report panel or Return to conversation to return. The user guide remains available in the top navigation.

The image shows a mobile application interface for a 'Laboratory conversation' workspace. At the top, there's a header with a 'Laboratory conversation' title and a 'Report' button. Below this is a disclaimer: 'Synthetic business data for demonstration only; no real service is provided.' The main area is titled 'What would you like to understand?' with a subtext 'Ask about a result, or attach your laboratory report.' There are two input fields: 'Understanding haemoglobin' and 'Creatinine and eGFR', both with a right arrow icon. Below these is a section titled 'Your result or question' with a text input area and a 'Send' button. At the bottom, there's an 'Attach report' button and a 'Send' button. A footer section contains instructions: 'Use the reference range printed on your report.', 'Examples use synthetic values. No diagnosis or treatment advice.', 'Educational information. No diagnosis, prescriptions or treatment changes.', and a link to 'View calculation rules'.

Laboratory conversation Report

Synthetic business data for demonstration only; no real service is provided.

What would you like to understand?

Ask about a result, or attach your laboratory report.

Understanding haemoglobin ↗

Creatinine and eGFR ↗

Your result or question

Type a question, or paste values with their units and reference ranges...

Attach report Send →

Use the reference range printed on your report. [Load example](#)

Examples use synthetic values. No diagnosis or treatment advice.

Educational information. No diagnosis, prescriptions or treatment changes.

[View calculation rules](#)

Actual application at a 390px viewport. Human usability testing has not yet been performed.

Review a report on a phone

The report panel fills the screen. Scroll inside it to inspect the image and all extracted fields, then select Confirm values. Keyboard focus stays inside the panel while it is open. Escape or Close report panel returns to the control that opened it. A different report requires a new analysis.

Documentation mock · Synthetic OCR values

Your report

Check the image and extracted values.

×

Upload a laboratory report

JPEG or PNG, up to 3 MB. Remove personal details first. PDF is not supported.

Choose File

No file chosen

📄

synthetic-report.png
34.9 KB · PNG

Review and edit the values

SYNTHETIC LABORATORY REPORT

Documentation fixture only. No patient information.

Test	Value	Unit	Reference range
Haemoglobin	10.8	g/dL	12 - 16
MCV	72	fL	80 - 100
Ferritin	7	ng/mL	15 - 150

Ranges are invented for demonstration, not clinical guidance.

Review values extracted from the image

Text from the image is unverified data. Edit unclear values before confirming.

Hb

Value

10.8

Unit

g/dL

Reference range

12-16

Actual mobile document panel with a MOCK OCR response. Values and report are synthetic.

Clear limits, throughout the workflow.

Laboratory information only. No diagnosis, prescriptions or treatment changes. Public reference ranges do not replace the ranges on your report.

Use de-identified or synthetic reports for demonstrations. When live access is enabled, relevant content may be sent to the configured provider. Ask the operator about the data policy before using real information.

PDF OCR, HIS and pharmacy integration, durable cloud settings and production release are not established by this guide.